

# Neop | 600 x 400 emergency trolleys with telescopic drawers

## Material and products immediately available for first aid

- Robustness: excellent resistance, one-piece structure made of aluminium
- 100 % hygiene: designed with a minimum of retention zones, and resistant to disinfecting products of hospital environment
- Customization thanks to a wide range of drawers and accessories
- Ward name identification
- Ergonomic and handy: designed for medical staff comfort
- Unique design and harmonious colours



Robust



Hygienic



Ergonomic



Customizable

### Colour

Garance red

- Height 100 cm
- Ergonomic gripping area
- Door wing lock with 1 batch of 20 seals. Access to the content of the trolley breaking the seal by simple twist. After each use, the trolley is reequipped, then closed by a single-use numbered seal.
- Drawers with full extension slides, with shock absorber and end stops, with perforated trays
- 5-drawer or 6-drawer version
- Frame for ISO trays adapted to most existing ISO trays
- Removable ISO tray for easy cleaning
- Cardiac massage board on support




## Standard configuration - PREMIUM

### EMERGENCY TROLLEY EQUIPPED WITH 5 telescopic drawers with compartmentalization

Treatment trolley 600 x 400 H100 equipped with :

- 1 hinged closure with seal securing
- 1 batch of 20 seals
- 1 telescopic drawer F600 H80 with full lip tray H50 equipped with a compartmentalization kit 16 boxes
- 3 telescopic drawers F600 H160 with openwork lip tray H100 equipped with a compartmentalization kit 9 boxes
- 1 telescopic drawer F600 H240 with openwork lip tray H200
- 1 cardiac massage board support
- 1 cardiac massage board
- 1 monitoring support
- 1 IV pole with 2 hooks
- 1 lateral standard rail – bottom right
- 1 oxygen tank holder

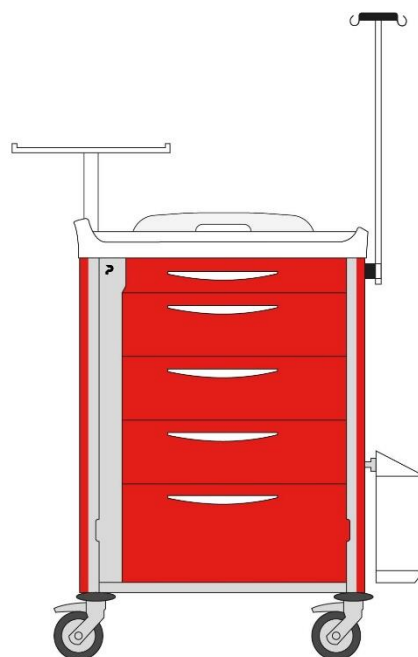


Reference	Description
8Y164H100T.U2P	Emergency trolley NEOP equipped 600x400 H100 - 5 telescopic drawers with compartmentalization

### EMERGENCY TROLLEY EQUIPPED WITH 5 telescopic drawers

Treatment trolley 600 x 400 H100 equipped with:

- 1 hinged closure with seal securing
- 1 batch of 20 seals
- 1 telescopic drawer F600 H80 with full lip tray H50
- 3 telescopic drawers F600 H160 with openwork lip tray H100
- 1 telescopic drawer F600 H240 with openwork lip tray H200
- 1 cardiac massage board support
- 1 cardiac massage board
- 1 monitoring support
- 1 IV pole with 2 hooks
- 1 oxygen tank holder
- 1 lateral standard rail – bottom right

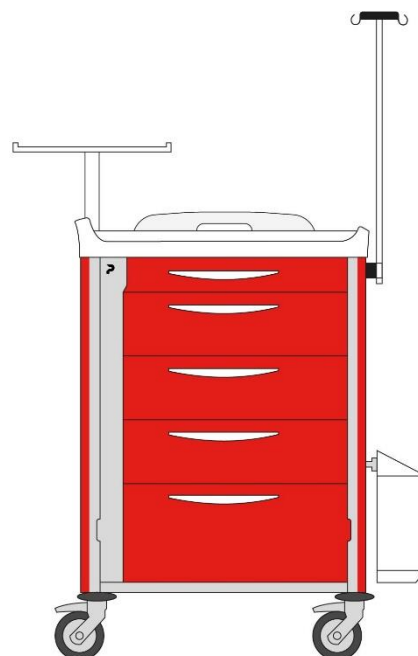


Reference	Description
8Y164H100T.U2	Emergency trolley NEOP equipped 600x400 H100 - 5 telescopic drawers

*Telescopic drawers fitted with lip trays can be divided with partitioning kits or combs, to be ordered separately.*

## EMERGENCY TROLLEY EQUIPPED WITH 6 telescopic drawers

Treatment trolley 600 x 400 H100 equipped with :

- 1 hinged closure with seal securing
- 1 batch of 20 seals
- 3 telescopic drawers F600 H80 with full lip tray H50
- 2 telescopic drawers F600 H160 with openwork lip tray H100
- 1 telescopic drawer F600 H240 with openwork lip tray H200
- 1 cardiac massage board support
- 1 cardiac massage board
- 1 monitoring support
- 1 IV pole with 2 hooks
- 1 lateral standard rail – bottom right
- 1 oxygen tank holder

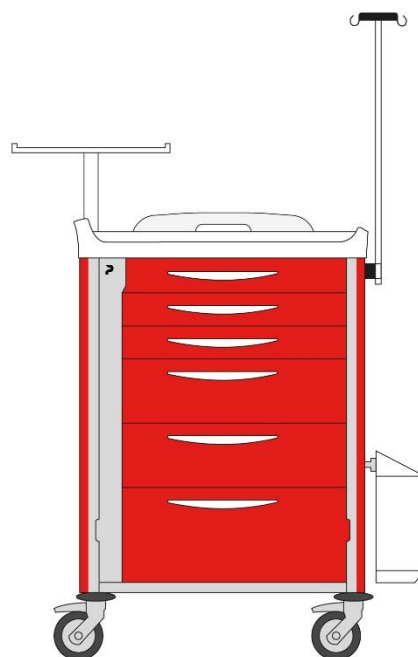


Reference	Description
8Y164H100T.U4	Emergency trolley NEOP equipped 600x400 H100 - 6 telescopic drawers

*Telescopic drawers fitted with lip trays can be divided with partitioning kits or combs, to be ordered separately.*

## BASIC EMERGENCY TROLLEY EQUIPPED WITH 5 telescopic drawers

Treatment trolley 600 x 400 H100 equipped with :

- 1 hinged closure with seal securing
- 1 batch of 20 seals
- 1 telescopic drawer F600 H80 with full lip tray H50
- 3 telescopic drawers F600 H160 with openwork lip tray H100
- 1 telescopic drawer F600 H240 with openwork lip tray H200
- 1 cardiac massage board support
- 1 cardiac massage board

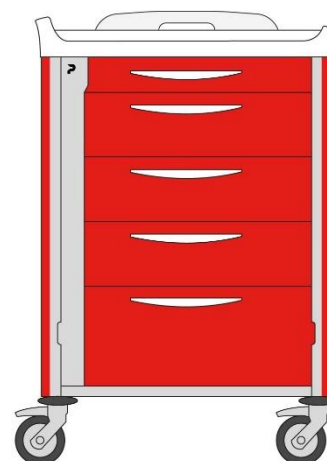


Reference	Description
8Y164H100T.U1	Basic emergency trolley NEOP 600x400 H100 - 5 telescopic drawers

*Telescopic drawers fitted with lip trays can be divided with partitioning kits or combs, to be ordered separately.*

## BASIC EMERGENCY TROLLEY EQUIPPED WITH 6 telescopic drawers

Treatment trolley 600 x 400 H100 equipped with :

- 1 hinged closure with seal securing
- 1 batch of 20 seals
- 3 tiroirs télescopiques F600 H80 with full lip tray H50
- 2 telescopic drawers F600 H160 with openwork lip tray H100
- 1 telescopic drawer F600 H240 with openwork lip tray H200
- 1 cardiac massage board support
- 1 cardiac massage board

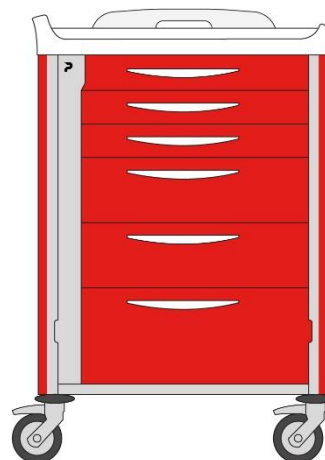

Reference	Description
8Y164H100T.U3	Basic emergency trolley NEOP 600x400 H100 - 6 telescopic drawers

*Telescopic drawers fitted with lip trays can be divided with partitioning kits or combs, to be ordered separately.*

# Technical data

## Overall dimensions of the trolley

The overall dimensions are measured with the wheels in the least bulky position.



## Useful dimensions of the worktop



## Worktop height

H. 1057 mm

## Materials

Structure	Complete welded and EPOXY anti-corrosion treated one-piece structure
Work surface	ABS / PMMA thermoplastic
Strip of colour	70 micron calendered PVC vinyl film with acrylic adhesive
Slides	Ballbearing slides made of galvanized steel with shock absorber
Closing system	Door wing: Aluminium with EPOXY painting
Wheels	Simple non marking rubber wheels - Polyamide cover (PA)
Bumpers	4 non marking round stop blocks made of EVA (Ethylene Vinyl Acetate)

## Additional information

Wheels	125 mm diameter – 2 breaking wheels + 2 swivelling wheels - wheel and wheel-pin double blocking	
Empty weight	38 kg	
Capacity	U1 / U2 : - 1 drawer HT 80 mm - 3 drawers HT 160 mm - 1 drawer HT 240 mm compact	U3 / U4 : - 3 drawers HT 80 mm - 2 drawers HT 160 mm - 1 drawer HT 240 mm compact
Locking system	Seals	
Colors	70 micron calendered polymeric PVC film coated with acrylic adhesive	
Equipment	Bare model (U1/U3) : Cardiac massage board	Model (U2/U4) with equipments: Cardiac massage board/oxygen tank holder / IV pole / monitoring support

## Packaging

Packing	Individually on pallet strapped and film-wrapped, or on wheels, with proper protection	Dimensions L X D x H	L 0,80 x P 1,20 x H 1,45 m
Gross Weight (trolley + packaging)	50 kg	Volume	1,15 m <sup>3</sup>